

## Volumetric Calculations

1.

15.5 mL of dilute hydrochloric acid reacts completely with 20.8 mL of  $0.0100 \text{ mol L}^{-1}$  calcium hydroxide. Calculate the concentration of the hydrochloric acid.

2.

20.5 mL of  $0.200 \text{ mol L}^{-1}$  hydrochloric acid reacts completely with 18.8 mL of sodium hydroxide solution. Calculate the concentration of the sodium hydroxide.

3. A  $0.100\text{ mol L}^{-1}$   $\text{H}_2\text{SO}_4$  solution is neutralised with  $10.0\text{ mL}$  of a solution of  $0.300\text{ M KOH}$ .
- a** Write a balanced equation for this reaction.
  - b** What volume of sulfuric acid was neutralised?

4.

15.0 mL of a nitric acid solution is required to react completely with 10.0 mL of a  $0.100 \text{ mol L}^{-1}$   $\text{Ca}(\text{OH})_2$  solution.

- a** Write a balanced equation for this reaction.
- b** What is the concentration of the nitric acid solution?

5.

18.26 mL of dilute nitric acid reacts completely with 20.00 mL of  $0.09927 \text{ mol L}^{-1}$  potassium hydroxide solution.

- a** Write a balanced chemical equation for the reaction between nitric acid and potassium hydroxide.
- b** Calculate the amount, in mol, of potassium hydroxide consumed in this reaction.
- c** What amount, in mol, of nitric acid reacted with the potassium hydroxide in this reaction?
- d** Calculate the concentration of the nitric acid.

6. Potassium hydrogen phthalate ( $\text{KH}(\text{C}_8\text{H}_4\text{O}_4)$ ) is used as a primary standard for the analysis of bases. Calculate the concentration of a standard solution prepared in a 50.00 mL volumetric flask by dissolving 2.042 g of potassium hydrogen phthalate in deionised water. The molar mass of  $\text{KH}(\text{C}_8\text{H}_4\text{O}_4)$  is  $204.2 \text{ g mol}^{-1}$ .

7.

Calculate the mass of anhydrous sodium carbonate ( $\text{Na}_2\text{CO}_3$ ) required to prepare 250.0 mL of a  $0.500 \text{ mol L}^{-1}$  standard solution.

8.

The concentration of a solution of barium hydroxide ( $\text{Ba}(\text{OH})_2$ ) was determined by titration with a standard solution of hydrochloric acid.

A 10.00 mL aliquot of  $\text{Ba}(\text{OH})_2$  solution was titrated with a  $0.125 \text{ mol L}^{-1}$  solution of HCl. Titres of 17.23 mL, 17.28 mL and 17.21 mL of HCl were required to reach the end point.

What is the concentration of the barium hydroxide solution?

9.

Potassium hydrogen phthalate,  $\text{KH}(\text{C}_8\text{H}_4\text{O}_4)$ , is a good primary standard for standardising alkali solutions. It contains one acidic hydrogen per formula unit. Potassium hydrogen phthalate (0.917 g) was dissolved in water and titrated with approximately  $0.2 \text{ mol L}^{-1}$  sodium hydroxide solution; 27.2 mL hydroxide solution was needed to reach the endpoint. Calculate the accurate molarity of the hydroxide solution.


**10.**

The acidity of a particular white wine was determined by titrating 25.0 mL of the wine with 0.0511 mol L<sup>-1</sup> sodium hydroxide solution; 8.70 mL was required. Calculate the molarity of hydrogen ions in the wine. Assume that the hydrogen ions come entirely from diprotic tartaric acid, H<sub>2</sub>C<sub>4</sub>H<sub>4</sub>O<sub>6</sub>, and calculate the concentration of tartaric acid in the wine in grams per 100 mL.



- a** Oxalic acid dihydrate,  $(\text{COOH})_2 \cdot 2\text{H}_2\text{O}$ , can be used as a primary standard for standardising alkali solutions. 0.291 g diprotic oxalic acid required 18.2 mL of a potassium hydroxide solution for exact neutralisation. Calculate the molarity of the hydroxide solution.
- b** 5.267 g anhydrous sodium carbonate was dissolved in water in a volumetric flask and the volume made up to 250 mL. 10 mL of this solution was pipetted into a conical flask and titrated with hydrochloric acid. 21.3 mL was needed to reach the equivalence point. Calculate the molarity of the hydrochloric acid solution. This solution was then used to determine the concentration of an unknown barium hydroxide solution. 25 mL of the barium hydroxide solution required 27.1 mL hydrochloric acid solution for exact neutralisation. Calculate the molarity of the barium hydroxide solution. In addition, calculate its concentration in grams per litre.

### Titration Involving Dilution

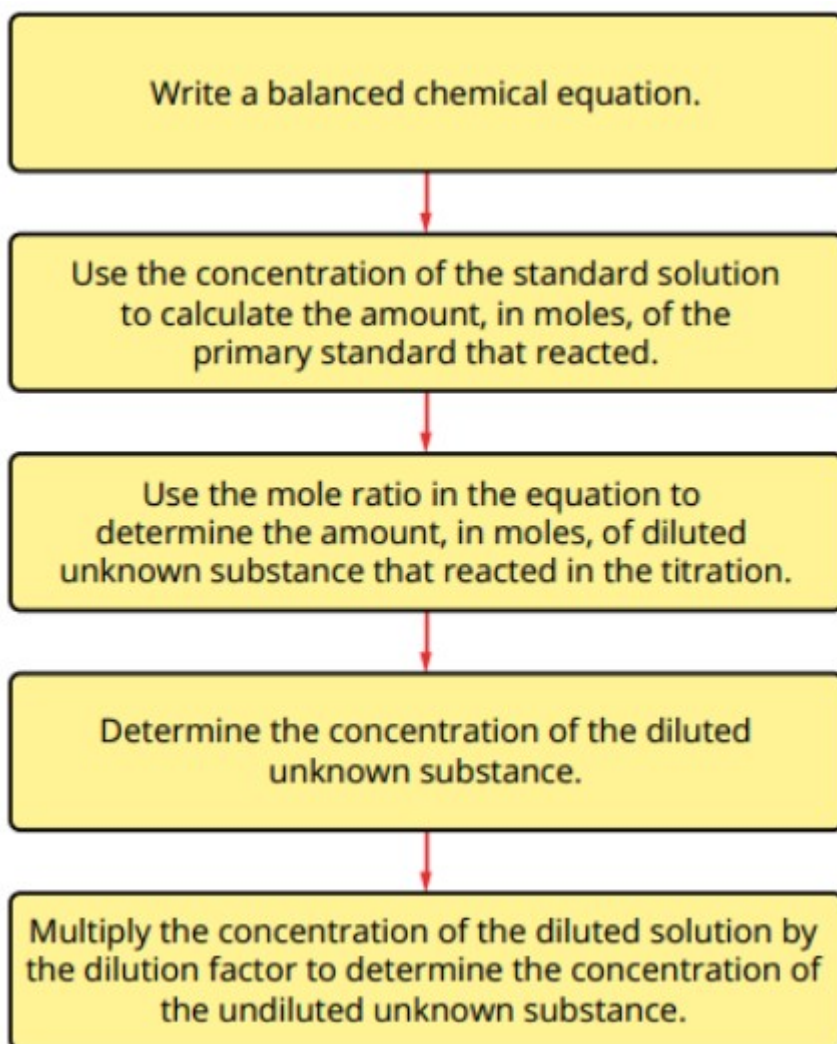

1.

A commercial concrete cleaner contains hydrochloric acid. A 25.00 mL volume of cleaner was diluted to 250.0 mL in a volumetric flask.

A 20.00 mL aliquot of 0.4480 mol L<sup>-1</sup> sodium carbonate solution was placed in a conical flask. Methyl orange indicator was added and the solution was titrated with the diluted cleaner. The indicator changed permanently from yellow to red when 19.84 mL of the cleaner was added.

Calculate the concentration of hydrochloric acid in the concrete cleaner.

2.

A commercial concrete cleaner contains hydrochloric acid. A 10.00 mL volume of cleaner was diluted to 250.0 mL in a volumetric flask.

A 20.00 mL aliquot of 0.2406 mol L<sup>-1</sup> sodium carbonate solution was placed in a conical flask. Methyl orange indicator was added and the solution was titrated with the diluted cleaner. The indicator changed permanently from yellow to pink when 18.68 mL of the cleaner was added.

Calculate the concentration of hydrochloric acid in the concrete cleaner.