


TOPIC 3:

TITRATIONS

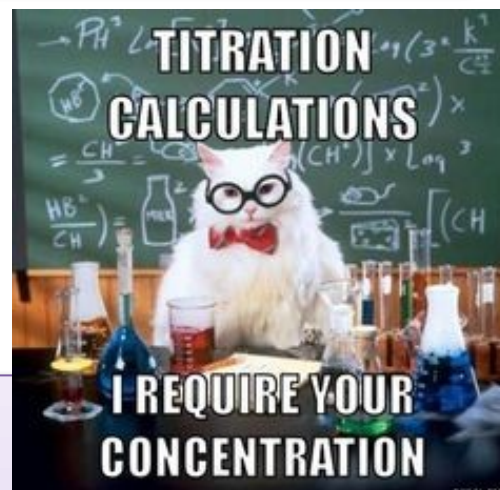

Syllabus dot-points

Science Understanding:

- ☐ volumetric analysis methods involving acid-base reactions rely on the identification of an equivalence point by measuring the associated change in pH, using appropriate acid-base indicators or pH meters, to reveal an observable end point
- ☐ data obtained from acid-base titrations can be used to calculate the masses of substances and concentrations and volumes of solutions involved

Science Investigation Skills:


- ☐ identify and distinguish between random and systematic errors, and estimate their effect on measured results

Timetable:

Wk	Monday P1	Tuesday P2	Wednesday P3	Thursday P5
1			Demo of NaOH vs HCl titration (pg 2-3) Standard solutions (pg 4-5) HW: Start calculations pg 24-39	Experiment: Making primary standard (pg. 6-7) Weekend HW: Finish calculations pg 24-39
2	Experiment: Standardising HCl solution (pg. 8-9) HW: Extended response preparation	Extended response task	Equivalence point and end point. Selecting appropriate indicators. (pg 10-13) HW: Read 'rinsing', 'errors', 'dilutions' (pg. 14-18)	Experiment: Diluting household ammonia and titrating against NH ₃ Weekend HW: Exploring Chem pg. 60-63 Start past exam questions
3	Chance to finish off calculations for NH ₃ experiment and seek help for difficult problems. Should be almost finished all HW and past exam Qs by this stage!!	Titration Practical Test (half of class)	Titration Practical Test (other half of class)	Titration Written Test

TOPIC 3:

TITRATIONS

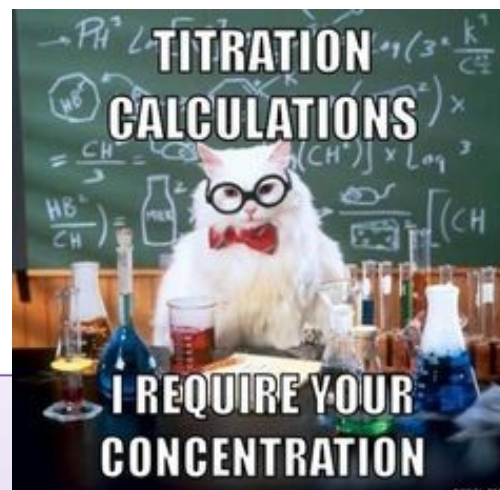

Syllabus dot-points

Science Understanding:

- ☐ volumetric analysis methods involving acid-base reactions rely on the identification of an equivalence point by measuring the associated change in pH, using appropriate acid-base indicators or pH meters, to reveal an observable end point
- ☐ data obtained from acid-base titrations can be used to calculate the masses of substances and concentrations and volumes of solutions involved

Science Investigation Skills:

- ☐ identify and distinguish between random and systematic errors, and estimate their effect on measured results

Timetable:

Wk	Tuesday P1	Wednesday P2	Thursday P3	Friday P5
1			Demo of NaOH vs HCl titration (pg 2-3) Standard solutions (pg 4-5) HW: Start calculations pg 24-39	Experiment: Making primary standard (pg. 6-7) Weekend HW: Continue calculations pg 24-39 Prepare for extended response
2	Extended Response Task	Experiment: Standardising HCl solution (pg. 8-9) HW: Finish calculations pg 24-39	Equivalence point and end point. Selecting appropriate indicators. (pg 10-13) HW: Read 'rinsing', 'errors', 'dilutions' (pg. 14-18)	Experiment: Diluting household ammonia and titrating against NH ₃ Weekend HW: Exploring Chem pg. 60-63 Start past exam questions
3	Chance to finish off calculations for NH ₃ experiment and seek help for difficult problems. Should be almost finished all HW and past exam Qs by this stage!!	Titration Practical Test	Titration Written Test	<i>Test results / Exam revision</i>


ACID-BASE TITRATIONS

Titration is a method for accurately determining the concentration of a substance.

There are two substances in a titration:

- A substance of known concentration (called a standard solution)
- A substance of unknown concentration

Today your teacher will demonstrate how this procedure could be used to find the concentration of a solution of HCl using a standard solution of NaOH .

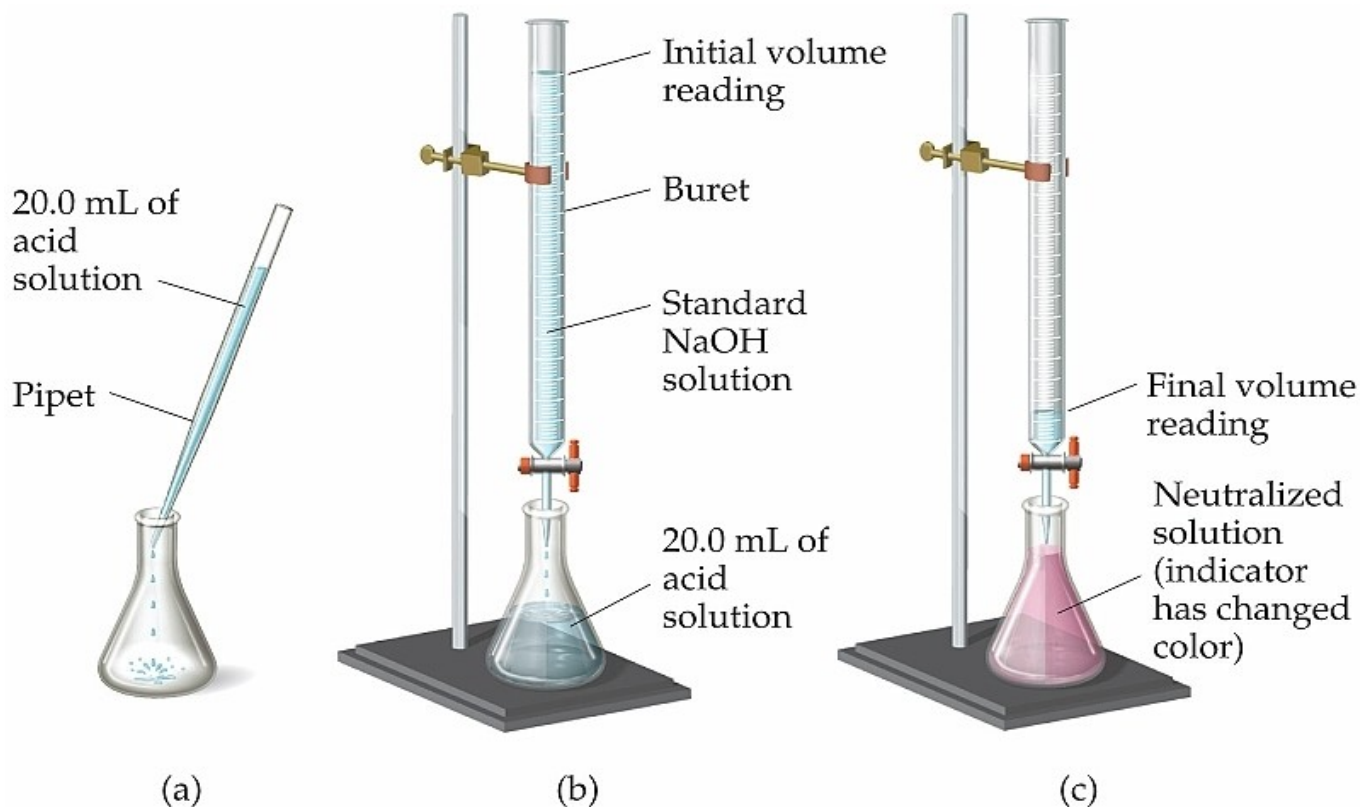

Copyright © 2009 Pearson Prentice Hall, Inc.

Step A:

- A volumetric pipette is rinsed with HCl .
- A conical flask is rinsed with distilled water.
- The pipette is used to draw up 20.0 mL of HCl , which is then added to the flask
- Two drops of an appropriate indicator is added

Step B:

- A burette is rinsed with NaOH solution
- NaOH is added to the burette until it is almost full
- Some NaOH is let out of the tap at the bottom to remove any air bubbles
- The initial volume of the burette is read
- The conical flask is placed under the burette

Step C:

- NaOH is added from the burette until the indicator has a permanent colour change
- The final volume of the NaOH solution is read

Measurements:

Titration is normally repeated until you have three measurements that are concordant ('concordant' means 'in agreement' or 'consistent'). This normally means three titre volumes that are within 0.1 mL of each other. Any inconsistent results are **not** included in the average titre volume.

	Trial 1	Trial 2	Trial 3	Trial 4
Initial volume (mL)				
Final volume (mL)				
Titre volume (mL)				

Average volume of NaOH: _____ mL

Concentration of NaOH solution: _____ mol L⁻¹ (from label on bottle)

- a) Calculate the moles of NaOH needed to neutralise the HCl. **Show working.**

.....

.....

.....

.....

.....

- b) The reaction occurring in the conical flask is: $\text{HCl} + \text{NaOH} \rightarrow \text{NaCl} + \text{H}_2\text{O}$
Calculate the moles of HCl in the conical flask. **Show working.**

.....

.....

.....

- c) Calculate the concentration of the original HCl solution. **Show working.**

.....

.....

.....


.....

.....

Standard Solutions

A **standard solution** is a solution with a **known concentration of substance**.


Standard solutions can be **primary standards** or **secondary standards**.

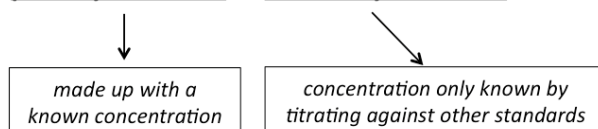

Primary Standards

A **primary standard** is substance that can be used to create a solution of known concentration.

This is performed by dissolving a known mass of the substance in a known volume of distilled water.

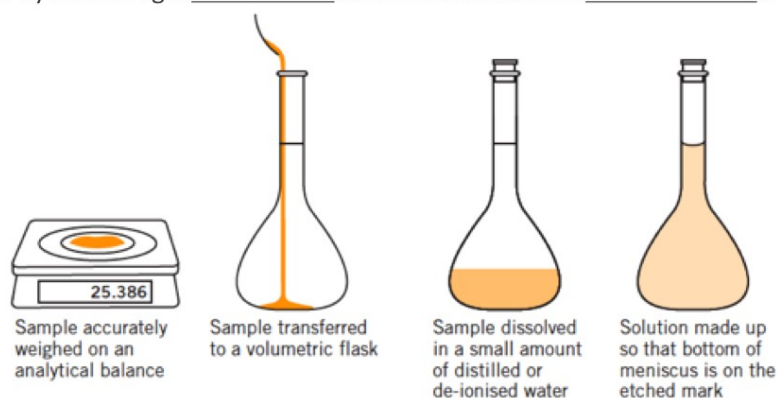

Example: Oxalic acid dihydrate ($\text{H}_2\text{C}_2\text{O}_4 \cdot 2\text{H}_2\text{O}$) can be used as a primary standard. As part of an experiment, 2.527 g of oxalic acid dihydrate was dissolved in 500 mL of water. Calculate the concentration of the oxalic acid solution.

$$M(\text{H}_2\text{C}_2\text{O}_4 \cdot 2\text{H}_2\text{O}) = 126.068 \text{ g mol}^{-1}$$

.....

.....

.....

.....

.....

.....

.....

.....

.....

Characteristics of Primary Standards

Primary standards have the following characteristics:

1. High purity
2. Stable in air (doesn't absorb or lose water moisture, doesn't react with gases in air)
3. Soluble in water
4. Has a relatively high molar mass

Appropriate primary standards:

- Anhydrous sodium carbonate (Na_2CO_3)
- Hydrated oxalic acid ($\text{H}_2\text{C}_2\text{O}_4 \cdot 2\text{H}_2\text{O}$)

Inappropriate primary standards:

Substance	Why it can't be used as a primary standard
Sodium hydroxide (NaOH)	Absorbs water moisture from air (hygroscopic) - Can absorb so much water from air that it will even dissolve (deliquescence) Reacts with CO_2 in air to form sodium carbonate
Concentrated HCl	Volatile. Loses HCl as gas fumes.
Concentrated H_2SO_4	Absorbs water moisture from air (hygroscopic)
<u>Hydrated</u> sodium carbonate	Loses water moisture to air (efflorescence)

The above issues mean you could not use **mass** to prepare a known **concentration** of solution.

Secondary Standards

Secondary standards are solutions whose concentration is determined by titration against other standards.

e.g. You want a standard solution of NaOH(aq) , but cannot prepare it as a primary standard.

Instead, prepare a solution of NaOH that is **approximately** 0.1 mol L^{-1} , and then titrate it against a primary standard solution of oxalic acid. This will tell you the **actual** concentration of the NaOH(aq) .

Preparing a Standard Soln of Na_2CO_3


On the following page are instructions for preparing a primary standard using Na_2CO_3 .

Before doing the experiment, calculate the mass of anhydrous Na_2CO_3 that would be needed to make up 500 mL of 0.05 mol L^{-1} solution. This will be the approximate mass you will try to measure out during the experiment.

Experiment: Preparation of a standard sodium carbonate solution

Background:

Standard solutions are solutions whose concentrations are accurately known. Volumetric techniques rely on accurate measurement of quantities. Mass and volume are the usual quantities measured. This experiment is the first in a series that prepares some standard solutions.

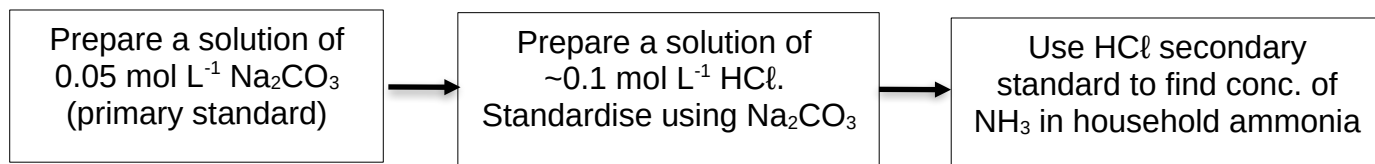

In this experiment you will prepare 500 mL of approximately 0.05 mol L⁻¹ Na₂CO₃ solution whose concentration is accurately known.

Equipment:

- Mass balance
- Volumetric flask (500 mL)
- Oven
- Desiccator
- Beaker (250 mL)
- Washbottle
- Storage bottle (approximately 500 mL)
- Distilled water
- Anhydrous sodium carbonate (Na₂CO₃) (~4 g)
- Stirring rod

Procedure:

1. **(Before the experiment)** Calculate the mass of anhydrous Na₂CO₃ required to make up 500 mL of 0.05 mol L⁻¹ solution.

Answer: _____ g


2. **(Before the experiment)** Your teacher will place the Na₂CO₃ in an oven at 270 °C for 30 minutes to remove any water, and then leave the anhydrous Na₂CO₃ in a desiccator to cool.
3. Accurately weigh out into a 250 mL beaker a mass of Na₂CO₃ approximately equal to that calculated. You should not waste time trying to weigh out exactly the mass calculated. Any value within ±0.5 g of the target will be fine, so long as you correctly record the actual mass used.

Mass of Na₂CO₃ used: _____ g

4. Dissolve the solid in ~100 mL of distilled water
5. Transfer this solution to a 500 mL volumetric flask. To ensure that all of the Na₂CO₃ has been transferred, rinse the beaker several times with about 20 mL portions of distilled water, adding each washing to the volumetric flask.

6. Make up the solution to precisely 500.0 mL with distilled water, until the meniscus of the water level is a few cm below the etched line on the flask.

Add the last few drops using a plastic pipette to ensure you do not overfill the flask. Stop when the meniscus is touching the etched line.

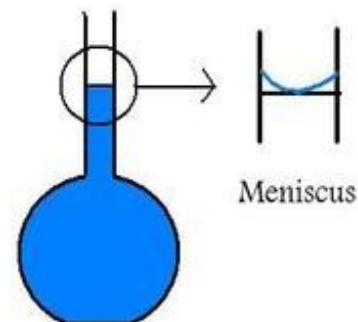

7. Place the stopper on the flask and invert the flask a few times to thoroughly mix the solution
8. Calculate the concentration of your standard solution:

.....


.....

.....

.....

.....

9. Use a small amount of the solution to rinse a storage bottle. Then, transfer your solution to the storage bottle. Label it with your name, the type of solution and its concentration.

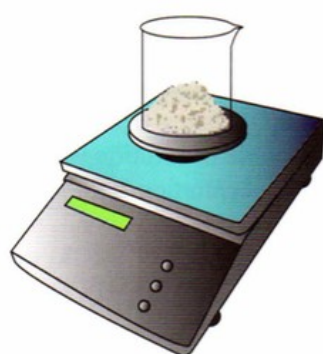

1. Weigh out sodium carbonate



2. Dissolve sodium carbonate in water



3. Transfer sodium carbonate solution to volumetric flask



4. Add water to make the solution up to the calibrated mark



5. Thoroughly mix the solution

Experiment: Standardisation of HCl solution

Background:

Hydrochloric acid cannot be used as a primary standard. In this experiment you will standardise an approximately 0.1 mol L^{-1} solution to determine its exact concentration by titration against the standard Na_2CO_3 solution prepared in the previous experiment. The equation for the reaction is:



Note that one mole of Na_2CO_3 reacts with two moles of HCl.


The carbon dioxide produced in this reaction results in the formation of carbonic acid (H_2CO_3) and therefore the reaction will have a somewhat acidic pH at the end of the reaction (pH ~3.5). Therefore an indicator that changes colour in this region must be used, such as methyl orange (yellow to red when pH changes from 4.4 to 3.1) or bromophenol blue (blue to yellow for pH change from 4.6 to 3.0).

Equipment:

- beakers (2x 100 mL, 1x 50 mL)
- burette and stand
- funnel
- pipette (20 mL)
- pipette filler
- conical flasks (4x 250 mL)
- storage bottle containing approximately 0.1 mol L^{-1} HCl
- storage bottle containing Na_2CO_3 primary standard from previous experiment
- methyl orange or bromophenol blue indicator (a few drops)
- distilled water (600 mL)

Procedure:

1. Label the 100 mL beakers as " Na_2CO_3 " and "HCl".
Label the 50 mL beaker as "waste"
2. Place about 100 mL of the standard Na_2CO_3 solution into the correctly labelled beaker. If the beaker is wet, rinse with a little of the Na_2CO_3 solution first.
3. Rinse the 'Rough' 250 mL conical flask with distilled water, and then rinse a clean 20 mL pipette with some of the Na_2CO_3 solution. Pipette a 20 mL aliquot of the Na_2CO_3 solution into the conical flask. Add 2-3 drops of indicator into the flask.
4. Place about 100 mL of HCl into a clean beaker. Again, if necessary, rinse the beaker with a little of the HCl solution first.
5. Rinse a clean burette with some of the HCl and then fill the burette with the solution. Pour some HCl into the waste beaker and ensure there are no bubbles in the bottom of the burette.
6. Note and record the level of acid in the burette. Obtain a rough estimate of the titration volume by running acid quickly from the burette while constantly swirling the liquid in the conical flask. Stop delivery of the acid as soon as a permanent colour change is obtained. Note and record the acid level in the burette and determine the approximate volume of acid required.

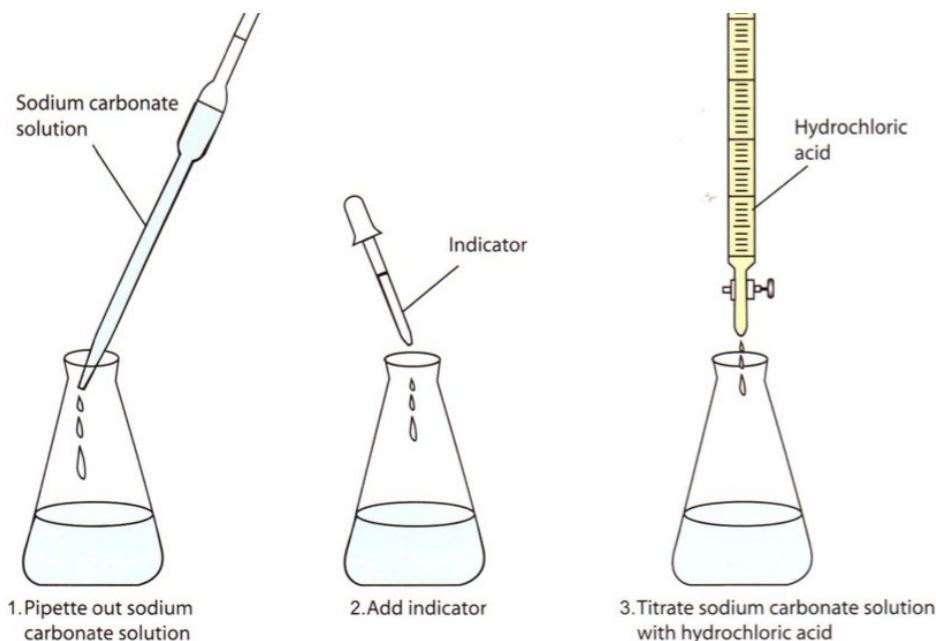

7. Prepare another conical flask containing 20 mL of the Na_2CO_3 solution and 2-3 drops of indicator. This time add the acid quickly from the burette with constant swirling of the flask, until the volume added is within 2-3 mL of the approximate volume required. Rinse the inside of the conical flask with a jet of water from a wash bottle to return any splashed solution to the bulk. Continue adding acid drop by drop, and with constant swirling, until the addition of one drop is sufficient to produce a permanent colour change. Note and record the level of the acid in the burette at the end point.
8. Repeat the accurate titration with further 20 mL portions of Na_2CO_3 solution until consistent titration volumes are obtained. These should be within 0.2 mL of each other.

Results

	Rough estimate	Accurate titrations			
		1	2	3	4 (if needed)
Initial volume (mL)					
Final volume (mL)					
Titre volume (mL)					

Calculate the concentration of the HCl solution:

.....

.....

.....

.....

.....

.....

.....

.....

.....

Indicators in Acid-Base Titrations

In an acid-base titration, an indicator needs to be selected that will change colour at the same time that the neutralisation reaction is finished.

Equivalence point: The point in a titration at which chemically equivalent amounts of acid and base have been added.

- For a 1:1 mole ratio (e.g. $\text{HCl} + \text{NaOH}$) this is when equal moles of acid and base have been added.
- For a 2:1 mole ratio (e.g. $\text{H}_2\text{SO}_4 + \text{NaOH}$) this is when $n(\text{H}_2\text{SO}_4) \text{ added} = 2 \times n(\text{NaOH})$

End point: The point in a titration at which the indicator changes colour

If an appropriate indicator has been picked, the equivalence point and the end point will be the same.

Common indicators.

Bolded = Remember!

Non-bolded = Here for reference only

Indicator	Colour on acidic side	Range of colour change	Colour on basic side
Methyl violet	Yellow	0.0 - 1.6	Violet
Bromophenol blue	Yellow	3.0 - 4.6	Blue
Methyl orange	Red	3.1 - 4.4	Yellow
Methyl red	Red	4.4 - 6.2	Yellow
Litmus	Red	5.0 - 8.0	Blue
Bromothymol blue	Yellow	6.0 - 7.6	Blue
Phenolphthalein	Colourless	8.3 - 10.0	Pink
Alizarin yellow	Yellow	10.1 - 12.0	Red


Calculate the change in pH when 0.100 mol L⁻¹ NaOH is added to 20.0 mL of 0.100 mol L⁻¹ NaOH

Important formulae: $n=cV$ $c=n/V$ $[\text{H}^+]=K_w/[\text{OH}^-]$ $[\text{OH}^-]=K_w/[\text{H}^+]$ $\text{pH}=-\log[\text{H}^+]$

V(HCl) initial (mL)	V(NaOH) added (mL)	Total volume (mL)	n(HCl) initial (mol)	n(NaOH) added (mol)	n(H ⁺) excess (mol)	n(OH ⁻) excess (mol)	[H ⁺] (mol/L)	[OH ⁻] (mol/L)	pH
20.0	0.0	20.0	0.00200	0.00000	0.00200	—	0.100	1.00×10^{-13}	1.00
20.0	5.0	25.0	0.00200			—			
20.0	10.0	30.0	0.00200			—			
20.0	15.0	35.0	0.00200			—			
20.0	19.0	39.0	0.00200			—			
20.0	19.5	39.5	0.00200			—			
20.0	19.8	39.8	0.00200			—			
20.0	20.0	40.0	0.00200	0.00200	—	—	1.00×10^{-7}	1.00×10^{-7}	7.00
20.0	20.2	40.2	0.00200		—				
20.0	20.5	40.5	0.00200		—				
20.0	21.0	41.0	0.00200		—				
20.0	25.0	45.0	0.00200		—				
20.0	30.0	50.0	0.00200	0.00300	—	0.00100	5.00×10^{-13}	0.0200	12.30

Graph:


Change in pH as 0.10 mol L^{-1} NaOH is added to 20.0 mL of 0.10 mol L^{-1} HCl

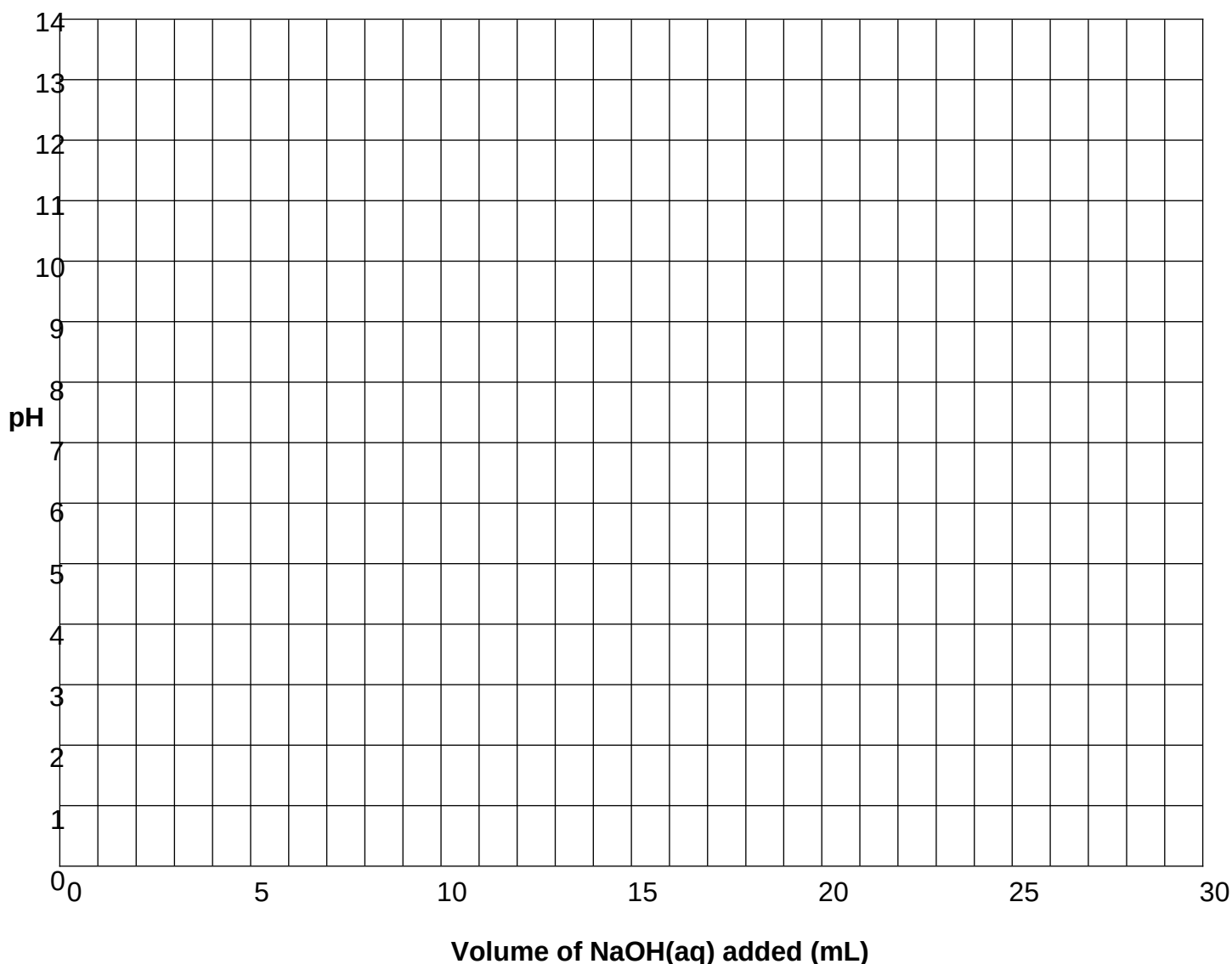

Questions.

1. On the graph, label the **equivalence point**.
2. What is the volume of NaOH added at the equivalence point? _____ mL
3. If you were using **phenolphthalein indicator**, at what **volume** would the indicator change colour from colourless to purple?
4. If you were using **methyl orange indicator**, at what **volume** would the indicator change colour from red to yellow?
5. On the basis of your answers above, explain whether phenolphthalein and/or methyl orange indicator would be suitable for this titration.

Weak Acid + Strong Base

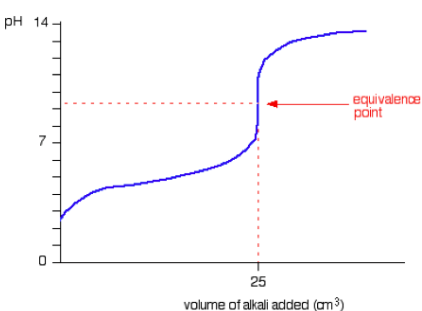

e.g. CH_3COOH and NaOH


- $$\text{CH}_3\text{COOH} + \text{NaOH} \rightarrow \dots\dots\dots$$

- Hint: This relates to something we learnt in Term 1 during the Acids & Bases topic...*

[illegible]

Choosing an Appropriate Indicator

Strong Acid + Weak Base e.g. HCl and NH_3

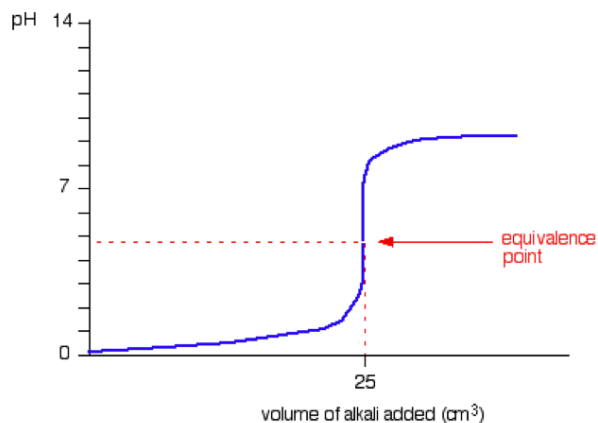


Which indicator would you use for this reaction:

phenolphthalein or methyl orange

Explain.

Weak Acid + Strong Base e.g. CH_3COOH + NaOH

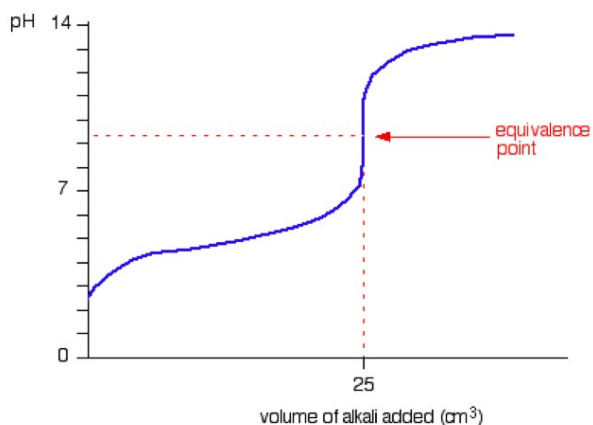


Which indicator would you use for this reaction:

phenolphthalein or methyl orange

Explain.

Weak Acid + Weak Base e.g. CH_3COOH + NH_3

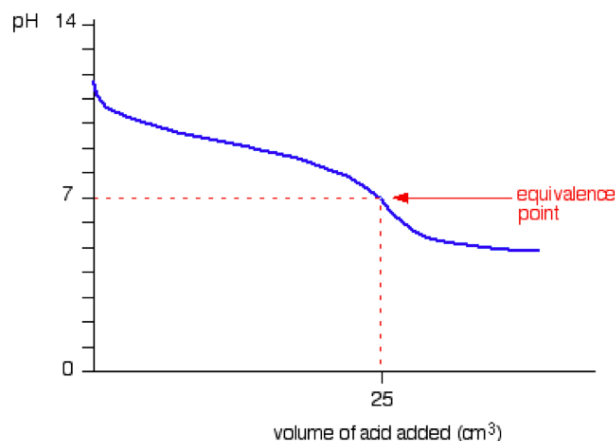

You will probably not come across a weak acid-weak base titration.

There is no 'steep part' on the titration curve for this combination. This makes it difficult to find the equivalence point using an indicator.

For these types of titrations you would normally use a pH meter and generate a pH graph like to one to the left to find the equivalence point.

SUMMARY:

Combination	Equivalence point	Appropriate indicator(s)
Strong acid + Strong base	Neutral	
Strong acid + Weak base	Acidic	
Weak acid + Strong base	Basic	

Rinsing

During the previous experiment, you had to rinse certain glassware with certain solutions.

Step	Reason
The pipette used to collect Na_2CO_3 was rinsed with a solution of Na_2CO_3	<p>If we had rinsed the pipette with water instead, then there would be some drops of water left in the pipette. This would have <u>diluted</u> the Na_2CO_3 solution.</p> <p>Part of the calculation to find $n(\text{HCl})$ was calculating: $n(\text{Na}_2\text{CO}_3) = c \times V$</p> <p>If the Na_2CO_3 had been diluted by water in the pipette then...</p> <ul style="list-style-type: none">• There would be less moles of Na_2CO_3 in the flask• Less volume of HCl would need to be added to the flask• We would <u>overestimate</u> the concentration of HCl
The conical flask was rinsed with distilled water before the 20.00 mL of Na_2CO_3 was added	<p>The purpose of the pipette was to collect <u>exactly</u> 20.00 mL of Na_2CO_3. We already know how many moles of Na_2CO_3 are in the pipette ($n = c \times V$).</p> <p>If we had rinsed the conical flask with Na_2CO_3 then...</p> <ul style="list-style-type: none">• There would be excess moles of Na_2CO_3 in the flask• More volume of HCl would need to be added to the flask• We would <u>underestimate</u> the concentration of HCl <p>Adding extra water to the conical flask would not in any way affect the number of moles of Na_2CO_3, and would have no effect on the results.</p>
The burette was rinsed with HCl	<p>Similar to the pipette, if we had rinsed the burette with water then there would be some drops of water left in the burette. This would have <u>diluted</u> the HCl solution.</p> <p>If the HCl solution was diluted then we would have needed a greater volume of HCl to react with the Na_2CO_3 solution. This means that we would have <u>underestimated</u> the concentration of HCl when calculating $c(\text{HCl}) = n / V$.</p>

As a general guideline:

- Pipettes, burettes and storage bottles are rinsed with the solution they will contain
- Conical flasks and volumetric flasks get rinsed with distilled water

It is important that you know what to rinse each piece of glassware with. One type of exam question might be to write a procedure for performing a titration. Your procedure should mention what each piece of glassware should be rinsed with.

Another type of exam question could be to explain the effects of incorrect rinsing procedures. Make sure that you understand the logic used in the above three examples, as this is the sort of mental working (and explanation) that you would need to be able to re-create in an exam.


Errors

In science, the word **error** does not mean the same thing as **mistake**. Experimental error means the difference between a **measurement** and the **true value**. All measurements have some form of error.

Accuracy and precision

- **Accuracy** refers to how close your value is to the true value.
- **Precision** refers to how consistent your results are.

The following image shows a depiction of accuracy and precision in archery.

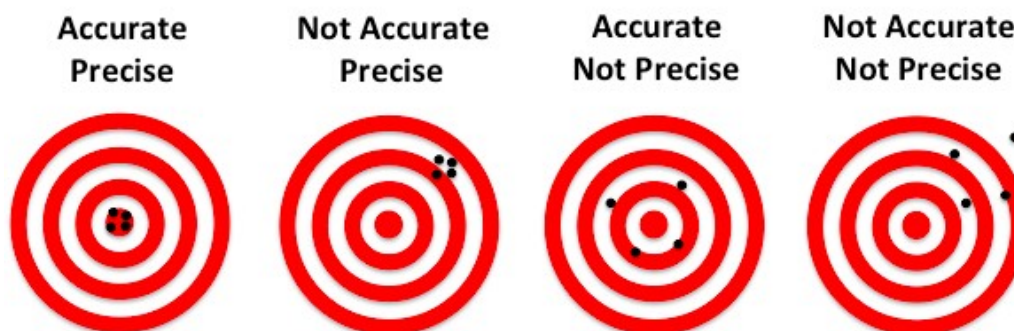

These concepts can be applied to titrations.

Imagine that for a particular titration, the correct titre volume is 18.21 mL.

The following are experimental results from three students:

	Titre Volume (mL)			
	Trial 1	Trial 2	Trial 3	Average
Student A	17.64	17.67	17.66	17.657
Student B	17.04	18.20	17.33	17.523
Student C	18.23	18.21	18.23	18.223

How would you describe the **accuracy** and **precision** of the three students?

Student A is:

Student B is:

Student C is:

Systematic errors

Systematic errors cause your results to be **consistently higher** or **consistently lower** than the correct value.

Examples of systematic errors might be:

- Always using the wrong indicator
- Always rinsing the pipette or burette with water
- Always measuring from the top of the meniscus when using a pipette


Systematic errors could explain how Student A achieved results that were precise but inaccurate. e.g. perhaps every time they rinsed the pipette with water before using it, which diluted the solution going into the conical flask and meant they needed less volume of solution from the burette.

Most systematic errors in titrations can be eliminated with proper technique.

Random errors

Random errors cause **inconsistent results**. Sometimes it will make your answer **higher** than the correct value, and sometimes it will make your result **lower** than the correct value.

Some random errors are caused by user error. For example, volumes of solution should be read at eye-level, but someone might sometimes look from above eye-level and look from below eye-level. This means sometimes the measured volume would be too large and other times it would be too small.

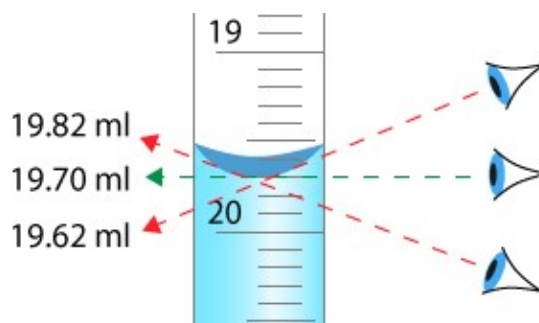

Random errors can never be fully eliminated, even when using information with the correct technique. There are limitations to the accuracy of equipment itself and the ability of the people. Consider measuring burette volumes – there is a limit to how accurately you can read volumes using the scale on the burette.

****IF**** all of the following equipment is used correctly...

- Burettes are able to measure volumes to the nearest 0.05 mL
- 250 mL volumetric flasks have a volume of 250.00 ± 0.12 mL
- 20 mL pipettes have a volume of 20.00 ± 0.06 mL
- A two decimal point mass balance can measure masses to nearest 0.01 g

The effects of random errors can be minimised by performing repeat trials, removing obvious outliers, and averaging the rest of the results.

Dilutions

Sometimes one of the solutions in a titration will be too concentrated, and will need diluting. In such cases you will often need to calculate the concentration of the original (undiluted) solution.

You can use the following formula to do this:

$$c_1 \times V_1 = c_2 \times V_2$$

c_1 and V_1 are the concentration and volume of the original (undiluted) solution, and c_2 and V_2 are the concentration and volume of the diluted solution.

This can be rearranged to give:

$$c_1 = c_2 \times \frac{V_2}{V_1}$$

Example:


A student analysed the ethanoic acid concentration in a commercial brand of vinegar. To do this, a 25.0 mL sample of the vinegar was diluted to 500.0 mL in a volumetric flask. 20 mL aliquots of the diluted vinegar were analysed by titration and found to have an ethanoic acid concentration of 0.0429 mol L⁻¹. Determine the concentration of ethanoic acid in the original undiluted vinegar.

$$c_{\text{eth. original}} = c_{\text{eth. dilute}} \times \frac{V_{\text{eth. dilute}}}{V_{\text{eth. original}}} = 0.0429 \times \frac{500}{25} = 0.858 \text{ mol L}^{-1}$$

Experiment: Ammonia content in household ammonia

Background:

Household ammonia is a solution containing about 3% by mass of ammonia (about 2 mol L⁻¹). In this experiment you will prepare a diluted solution of household ammonia and titrate it against your standardised HCl solution.

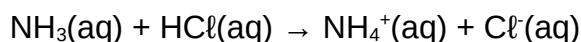

Equipment (Part A):

- Household ammonia solution
- 10 mL volumetric pipette
- 50 mL beaker
- Distilled water (300 mL)
- Mass balance

Equipment (Part B):

- Standardised HCl solution
- 20 mL volumetric pipette
- Burette and stand
- Funnel
- Indicator (methyl orange or phenolphthalein)
- 2x 100 mL beaker
- 1x 50 mL beaker
- 4x 250 mL conical flask

Procedure: Preparing the diluted NH₃ solution

1. Fill a beaker with 30 mL of household ammonia solution
2. Rinse a 10 mL pipette with the household ammonia solution
3. Rinse a 250 mL volumetric flask with distilled water.
4. Weigh the volumetric flask, then add 10.0 mL of household ammonia solution using the volumetric pipette. **Reweigh the volumetric flask before adding any water.**
(This information will be used later to calculate the % by mass of ammonia)

	Mass (g)
Empty vol. flask	
Vol. flask + 10 mL household ammonia	
10 mL household ammonia (calculated)	

5. Fill the remainder of the flask with distilled water up to the mark. Invert the flask to mix.
6. *If this is happening in a different lesson to the next part of the procedure, you will need to store your diluted ammonia solution in a storage bottle. Rinse the storage bottle with some of the diluted ammonia solution before filling the bottle.*

Procedure: Titrating the diluted ammonia solution vs. standardised HCl

Perform the titration using 20.0 mL aliquots of diluted ammonia solution, with your standardised HCl solution in the burette. Prior to beginning the experiment, consider what you should use to rinse each piece of glassware in the experiment. Check with your teacher before proceeding.

Glassware	Rinse with...
Burette	
20.0 mL pipette	
Conical flask	

This will be a titration between a strong acid and a weak base. Which indicator will you use in this reaction?

Perform a rough titration and then two to four additional titrations. Stop when you have three concordant titres.

Results:	Rough estimate	Accurate titrations			
		1	2	3	4 (if needed)
Initial volume (mL)					
Final volume (mL)					
Titre volume (mL)					

Calculations:

Calculate the concentration of NH_3 in the original (undiluted) ammonia solution (in mol L^{-1}).

.....

.....

.....

.....

.....

.....

.....

.....

.....

.....

.....

.....

During the preparation of the diluted solution you should have measured the mass of 10.0 g of the original household ammonia solution. Using this value and your answer above, calculate the percentage by mass of NH_3 in the household ammonia solution.

.....

.....

.....

.....

.....

.....

.....

.....

.....

.....

Uncertainty Calculations (advanced)

This information is beyond the scope of the Year 12 Chemistry syllabus, however is useful additional information for high achieving students and/or those intending to pursue a Science degree at university.

Uncertainty values are a measure of the random error inherent in equipment. For example, Class B pipettes can deliver a volume of 20.00 ± 0.06 mL. The " ± 0.06 mL" is the uncertainty of the pipette. Regardless of how correctly you use the pipette, this is the limitation to its accuracy.

Equipment	Uncertainty
10.00 mL pipette	± 0.01 mL (Class A pipette) ± 0.02 mL (Class B pipette)
20.00 mL pipette	± 0.03 mL (Class A pipette) ± 0.06 mL (Class B pipette)
Burette	± 0.10 mL for a titre volume (final - initial) (± 0.05 mL for each reading)
250 mL volumetric flask	± 0.12 mL
500 mL volumetric flask	± 0.20 mL
2 decimal point mass balance	± 0.01 g for measured mass (final - initial)
3 decimal point mass balance	± 0.001 g for measured mass (final - initial)

Percentage Uncertainty:

Titration involves the use of multiple measurements, each with its own uncertainty. To find out the overall uncertainty on the final calculation, we first need to convert each uncertainty to a percentage uncertainty.

$$\text{Percentage uncertainty} = \frac{\text{uncertainty}}{\text{reading}} \times 100$$

For example, the percentage uncertainty for a 20.0 mL pipette = $0.06/20 \times 100 = 0.3\%$

By summing the percentage uncertainty for each measurement, you get the overall percentage uncertainty for the calculated concentration.

EXAMPLE:

Joey performs a titration using NaOH and HCl. Before performing the experiment he dilutes the HCl by pipetting 20.00 mL into a 250 mL flask, and then makes the volume up to 250 mL. He adds this solution to the burette. He then pipettes 20.00 mL of NaOH into a conical flask, and performs the titration. It takes an average of 14.20 mL of HCl to neutralise the NaOH. From this, he calculates the concentration of HCl is $0.12572 \text{ mol L}^{-1}$.

What is the uncertainty in the concentration of HCl?

Working on following page...

Equipment	Measurement	Percentage uncertainty
Pipette HCl	20.00 ± 0.06 mL	0.30%
Vol. flask	250.00 ± 0.12 mL	0.048%
Pipette NaOH	20.00 ± 0.06 mL	0.30%
Titre volume	14.20 ± 0.10 mL	0.70%
Total:		1.348%

1.348% of 0.12572 is **0.002 mol L⁻¹**. (*The uncertainty only gets written to 1 significant figure*)

Therefore, Joey should report his concentration as **0.126 ± 0.02 mol L⁻¹**.

In other words, the concentration of the HCl is between 0.124 mol L⁻¹ and 0.126 mol L⁻¹.

Extension task: Calculate the uncertainty of your Na₂CO₃ concentration, HCl concentration and your %mass of NH₃ in the household ammonia.

Part 1: Sodium carbonate solution

	Measurement	Percentage uncertainty
Mass of Na ₂ CO ₃		
500 mL vol flask	500 ± 0.20	0.04%
Total:		

Concentration of sodium carbonate solution: _____

Part 2: Hydrochloric acid solution

	Measurement	Percentage uncertainty
Conc. of Na ₂ CO ₃		
20 mL pipette Na ₂ CO ₃		
Titre volume HCl		
Total:		

Concentration of HCl solution: _____

Part 3: Household ammonia solution

	Measurement	Percentage uncertainty
10.0 mL conc NH ₃		
250 mL vol flask		
20.0 mL dilute NH ₃		
Titre volume HCl		
Concentration of HCl		
Mass of 10 mL NH ₃ soln		
Total:		

Percentage by mass of NH₃ in household ammonia: _____

What was the biggest source of error in this experiment? How could it be minimised?

TITRATIONS PRAC TEST

The titrations topic will be assessed over two lessons:

- A one hour **practical test** that will assess your ability to correctly perform a titration. You will be marked based on your experimental results and how closely they match the expected results.
- A one hour **written test** that will assess your broader understanding of volumetric analysis, including both explanation-type questions and calculation-type questions. One of the calculation-type questions will be related to your practical test. The style of questions in this assessment will be similar to those in past exam papers.

Component	Weight of yearly mark
Part 1: Practical skills	2.0%
Part 2: Written questions	4.0%
Total	6.0%

Practical test details:

The practical test will be run on the date shown on the **front cover of this booklet**. For the larger class, the practical test will run over two lessons.

The practical test will require you to prepare a diluted vinegar solution and titrate it against a standardised solution of NaOH. You will be assessed on the quality and accuracy of your results. Your accuracy will be based on how closely your results align to mine (Mr McKenna's). I will perform the same titration with the same solutions. You will receive points based on how close your titre volumes are to my own titre volumes. If your titre volumes fall within the same uncertainty range as mine then you will receive full marks for accuracy. Marks will decrease as you become further away from my score.

Example:

Mr McKenna has a titre volume of 17.2 ± 0.2 mL

Student A has a titre volume of 17.1 ± 0.2 mL Full marks for Student A

Student B has a titre volume of 17.6 ± 0.2 mL Full marks for Student B
(as perhaps the 'true' value is 17.4 mL)

Student C has a titre volume of 16.5 ± 0.2 mL Not full marks for Student C

The only calculation you will need to do is to calculate your average titre volume. You will **not** need to perform any other calculations during the practical lesson, however there will be calculation questions on the written test related to this experiment.

The front page of the practical test is shown to the right.

PRACTICAL TEST – TITRATIONS

Section	Marks
Can fill volumetric flask up to the line	/ 1
Can consistently reach an end point	/ 2
Can accurately use a 20.0 mL pipette	/ 3
Accuracy and precision of titre volumes	/ 9
TOTAL:	/ 15

Overview:

In this practical test you will prepare a diluted solution of acetic acid from commercial vinegar. The acetic acid will be titrated against standardised sodium hydroxide.

Solutions supplied:

- A standardised solution of NaOH (150 mL)
- Commercial vinegar solution (50 mL - undiluted)
- Phenolphthalein indicator
- Distilled water

Supplied equipment:

- 250 mL volumetric flask
- 3x 100 mL beaker (for storing NaOH, conc. vinegar, dilute vinegar)
- 1x 50 mL beaker (waste)
- 20 mL volumetric pipette and filler
- Conical flasks
- Burette and stand
- Funnel

Part 1: Dilution of vinegar

Commercial vinegar is much more concentrated than the sodium hydroxide we will use, so it will require dilution before it can be titrated. Use the provided volumetric pipette and volumetric flask to prepare a diluted vinegar solution.

IMPORTANT: During this titration you will need to show the teacher a volumetric flask with 250.0 mL of solution. It is recommended that you do this while making up your solution during Part 1. Show the full volumetric flask before inverting the flask as some of the liquid will stick to the stopper during the inversion process.


Part 2: Titration against sodium hydroxide

This titration will be performed with the sodium hydroxide solution in the burette and 20 mL of diluted vinegar in a conical flask. Perform a number of titrations until you have achieved consistent results.

IMPORTANT: During this titration you will need to show the teacher two conical flasks which are at the correct endpoint for your chosen indicator. You should show both conical flasks at the same time so that the consistency of your endpoint can be evaluated.

Question 3:

A solution of NaOH was standardised by reacting with a $0.1008 \text{ mol L}^{-1}$ solution of H_2SO_4 . An average volume of 18.64 mL of sulfuric acid was required to neutralise the 20.00 mL aliquots of NaOH. Find the concentration of NaOH.

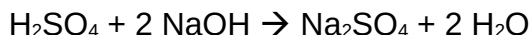


HINT: This problem is different from the others because sulfuric acid is a **diprotic acid**. There is a 2:1 mole ratio between H_2SO_4 and NaOH . $n(\text{NaOH}) = 2 \times n(\text{H}_2\text{SO}_4)$.

A: 0.1879 mol L⁻¹

Question 4:

A solution of phosphoric acid (H_3PO_4), a triprotic acid, was standardised by reacting with a $0.08555 \text{ mol L}^{-1}$ solution of NaOH . An average volume of 27.55 mL of phosphoric acid was required to neutralise the 20.00 mL aliquots of NaOH . Find the concentration of phosphoric acid.

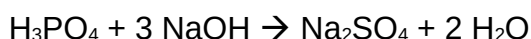

HINT: Once again, be careful of the mole ratio. It is not 1:1.

A: 0.02070 mol L⁻¹

Question 5:

The main ingredient in vinegar is acetic acid (CH_3COOH), a monoprotic acid. A solution of acetic acid was standardised by reacting with a $0.09782 \text{ mol L}^{-1}$ solution of NaOH . An average volume of 18.64 mL of NaOH was required to neutralise the 20.00 mL aliquots of CH_3COOH . Find the concentration of acetic acid.

HINT: No equation was written for this reaction. You should write a balanced reaction to find the mole ratio.

A: 0.09117 mol L⁻¹

Question 6:

The main ingredient in wine is tartaric acid ($\text{C}_4\text{H}_6\text{O}_6$), a diprotic acid. A solution of tartaric acid was standardised by reacting with a $0.1024 \text{ mol L}^{-1}$ solution of NaOH. An average volume of 18.64 mL of tartaric acid was required to neutralise the 20.00 mL aliquots of NaOH. Find the concentration of tartaric acid.

..... A: 0.05494 mol L⁻¹

Question 7:

A titration was performed using 20 mL aliquots of 0.04056 mol L⁻¹ solution of sodium carbonate. The sodium carbonate was titrated against a solution of nitric acid with a concentration of approximately 0.1 mol L⁻¹. The following information was recorded for the titration.

	Volume of HNO ₃ solution required			
	1	2	3	4
Initial volume (mL)	0.45	1.20	0.05	12.40
Final volume (mL)	20.80	20.40	19.25	31.65
Titre volume (mL)	20.35	19.20	19.20	19.25

Calculate the concentration of the nitric acid solution.

Hint: You will need to calculate the average titre volume to use in your calculation. Don't include any outliers in your average. Only include concordant (consistent) values. i.e. within ± 0.1 mL.

A: 0.08443 mol L⁻¹

Question 8:

A solution of sodium hydroxide was standardised by titrating against $1.067 \text{ mol L}^{-1} \text{ H}_2\text{SO}_4$. The titration was performed using 20 mL aliquots of H_2SO_4 with the NaOH in the burette. The following information was recorded for the titration:

	Volume of NaOH solution required				
	1	2	3	4	5
Initial volume (mL)	0.20	1.65	1.60	17.20	12.40
Final volume (mL)	17.85	18.05	17.90	36.15	28.70
Titre volume (mL)					

Calculate the concentration of the sodium hydroxide solution.

A: 2.613 mol L⁻¹

Question 9:

A solution of household ammonia was analysed by titration to determine the concentration of ammonia in the solution. Prior to analysis, the ammonia was **diluted** by pipetting a 10.00 mL aliquot into a 500 mL volumetric flask and making up to the mark using distilled water. The titration was then performed using 20.00 mL aliquots of the **diluted** ammonia solution in the conical flask, and 0.07844 mol L⁻¹ hydrochloric acid in the burette.

The following information was recorded for the titration:

	Volume of HCl solution required			
	1	2	3	4
Initial volume (mL)	0.75	7.20	4.30	4.95
Final volume (mL)	23.40	29.80	26.30	27.50
Titre volume (mL)				

Calculate the concentration (in mol L⁻¹) of ammonia in the **original (undiluted)** solution.

HINT: First, calculate the concentration of diluted ammonia using information about the titration. Then you can use the formula $c_1 \times V_1 = c_2 \times V_2$ to calculate the concentration of the undiluted ammonia solution.

..... A: 4.432 mol L⁻¹

A solution of vinegar was analysed by titration to determine the concentration of acetic acid. Prior to analysis, the vinegar was diluted by pipetting a 25.00 mL aliquot into a 250 mL volumetric flask and making up to the mark using distilled water. The titration was then performed using 20.00 mL aliquots of the 0.1005 mol L⁻¹ NaOH in the conical flask, and diluted vinegar in the burette.

	Volume of diluted CH ₃ COOH solution required			
	1	2	3	4
Initial volume (mL)	0.25	4.55	12.60	6.60
Final volume (mL)	20.30	24.25	32.30	26.20
Titre volume (mL)				

A: 1.022 mol L⁻¹


Question 11:

A solution of Na_2CO_3 was prepared by dissolving an unknown mass of anhydrous Na_2CO_3 into 250.0 mL of water. 20.00 mL aliquots of the resulting solution were titrated against $0.09835 \text{ mol L}^{-1} \text{ HCl}$. The average titre volume was 25.87 mL.

Calculate the mass of Na_2CO_3 dissolved in the volumetric flask.

HINT: Use your normal titration calculations to find the number of moles of Na_2CO_3 in the conical flask, and then from this find the number of moles of Na_2CO_3 in the volumetric flask.

Remember, a simple diagram of the steps can aid in understanding what is happening.

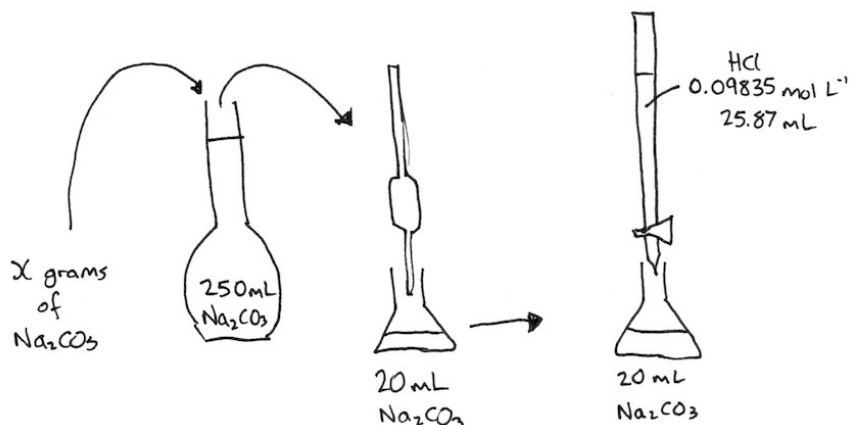


Question 12:

A solution of oxalic acid ($\text{H}_2\text{C}_2\text{O}_4$) was prepared by dissolving an unknown mass of oxalic acid dihydrate ($\text{H}_2\text{C}_2\text{O}_4 \cdot 2\text{H}_2\text{O}$) in a 500 mL volumetric flask. The oxalic acid was placed in the burette, and was titrated against 20.00 mL aliquots of $0.007930 \text{ mol L}^{-1}$ NaOH solution. The average titre volume was 22.40 mL.

Calculate the mass of oxalic acid dihydrate dissolved in the 500 mL volumetric flask.

Note: Oxalic acid ($\text{H}_2\text{C}_2\text{O}_4$) is a diprotic acid.

HINT: Use your normal titration calculations to find the number of moles of $\text{H}_2\text{C}_2\text{O}_4$ in the 22.40 mL titre volume, and from this find the number of moles of $\text{H}_2\text{C}_2\text{O}_4$ in the volumetric flask.

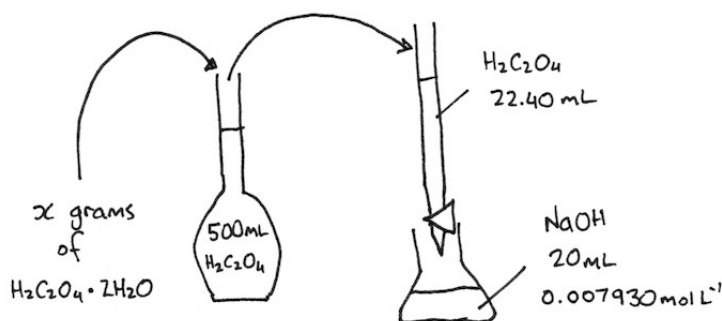

A: 0.2232 g

Question 13:

Washing soda is an impure mixture which mainly contains sodium carbonate decahydrate ($\text{Na}_2\text{CO}_3 \cdot 10\text{H}_2\text{O}$).

4.81 g of washing soda was dissolved in 500 mL of water and the flask was inverted multiple times to ensure thorough mixing. 20.00 mL aliquots of this solution were then titrated against 0.06544 mol L⁻¹ hydrochloric acid using methyl orange indicator. An average titre volume of 18.45 mL was required to reach the end point.

Calculate the percentage by mass of sodium carbonate decahydrate in the washing soda sample.

HINT: You will need to find the mass of $\text{Na}_2\text{CO}_3 \cdot 10\text{H}_2\text{O}$ in the sample using information from the titration, and then use the following formula:

$$\% Na_2CO_3 \cdot 10H_2O = \frac{m(Na_2CO_3 \cdot 10H_2O)}{m(\text{washing soda})} \times 100$$

A: 89.8%

Question 14:

Potassium hydroxide (KOH) cannot be used as a primary standard because it is deliquescent (absorbs water moisture from the air).

A solution of potassium hydroxide was prepared by dissolving 1.32 g of potassium hydroxide pellets in a 250 mL volumetric flask. A small amount of the solution was used to rinse a burette, and then the burette was filled with the potassium hydroxide solution.

Separately, 20.00 mL aliquots of 0.09667 mol L⁻¹ HCl were pipetted into clean conical flasks, and 2-3 drops of bromothymol blue indicator was added. KOH solution was added until the bromothymol blue changed colour from yellow to green. On average, 20.97 mL of KOH solution was required to reach the end point.

Calculate the %KOH and %H₂O in the potassium hydroxide pellets used to make the KOH solution.

A: 98.0%, 2.0%

The acidity of wine is due mainly to potassium tartrate (cream of tartar), a weak monoprotic acid with a molar mass of 188 g mol^{-1} . Three 50 mL samples of wine were titrated with 0.012 mol L^{-1} NaOH. The results of the titrations are shown below:

- Calculate the moles of potassium tartrate in 50 mL of wine
- Calculate the concentration of potassium tartrate in g L^{-1} .

(a) 0.0228 g; (b) 0.456 g L⁻¹

A solution of household vinegar was analysed to determine the acetic acid content. 10.00 mL of vinegar was added to a 250 mL volumetric flask and distilled water was added until the total volume of the solution was 250 mL. 20.00 mL aliquots of the diluted vinegar solution were pipetted into clean conical flasks, and 2-3 drops of phenolphthalein was added.

- Calculate the mass of CH_3COOH in 10 mL of the original undiluted vinegar
- Given that the density of household vinegar is 1.05 g/mL, calculate the percentage by mass of acetic acid in the household vinegar.

$$\%CH_3COOH = \frac{m(CH_3COOH)}{m(vinegar)} \times 100$$
$$\text{density} = \frac{\text{mass}}{\text{volume}}$$
[illegible]

A student was asked to determine the mass, in grams, of calcium carbonate present in a 0.125 g sample of chalk. The student placed the chalk in a 250 mL conical flask and added 50.00 mL of 0.200 mol L⁻¹ HCl using a volumetric pipette. The solution was stirred until no further reaction occurred.

Calculate the mass of calcium carbonate, in grams, present in the chalk sample.

$$n(\text{HCl reacted}) = n(\text{HCl initial}) - n(\text{HCl excess})$$

Find out how many moles of HCl reacted with the calcium carbonate, and use that to figure out how much calcium carbonate was in the chalk.

[illegible]

40

Question 18:

US federal regulations set the upper limit of ammonia (NH_3) in air in a work environment of $38.5 \mu\text{g}$ ($38.5 \times 10^{-6} \text{ g}$) per litre of air. To test the air quality in a manufacturing plant, 100.0 L of air was bubbled through 100 mL of $0.0105 \text{ mol L}^{-1} \text{ HCl}$.

10.0 mL aliquots of the reacted mixture is then titrated with $0.00588 \text{ mol L}^{-1}$ NaOH. On average, 13.10 mL of NaOH was required to neutralise the excess HCl.

Calculate the concentration of ammonia in the air sample (in $\mu\text{g/L}$) and state whether or not the manufacturer is in compliance with regulations.

HINT: This is another back titration. Figure out how many moles of HCl must have reacted with the NH_3 , and use that to find out how much NH_3 was in 100.0 L of air.

A: 48.6 $\mu\text{g/L}$, no

HOMEWORK:

Source	Set/Chapter	Questions to complete
This workbook	Calculation questions (pg. 24-39)	<div>1 2 3 4 5 6 7 8 9 10</div> <div>11 12 13 14 15 16 17 18</div>
Essential Chemistry	Set #7 (pg. 60-63)	<div>1 2 3 4 5 6 7 8 9 10</div> <div>11 12 13 14 15 16 17 18</div>
Past exam questions		<div>Page... 2 3 4 5 6 7 8 9</div> <div>10 11 12 13 14 15 16 17</div> <div>18 19 20 21 22 23 24 25</div> <div>26 27 28</div>